
Pulsebit EX

Personal ECG Tracker

Quick Guide

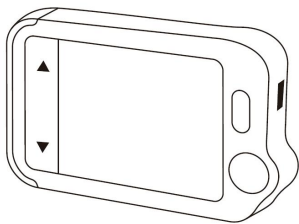Model: PB-20

English | Deutsche | Italiano | Español | Français

Table of contents

1. Disclaimer.....	2
2. Caution.....	1
3. Unpacking.....	1
4. About the Main Unit.....	2
5. ViHealth App.....	2
6. Desktop Software.....	3
7. Set Up.....	4
7.1 Charging.....	4
7.2 Power On/Off.....	4
7.3 Initial Settings.....	4
7.4 First Check.....	5
7.5 View Results.....	6
7.6 Reviewing Records.....	6
7.7 Syncing data to ViHealth App.....	7
7.8 Using the Software “Pulsebit Browser Pro”.....	10
8. Exploring More Features.....	11
8.1 Choosing Different ECG Types.....	11
8.2 Using the ECG Cable for Heart Checks.....	12
8.3 About Reusable Electrodes.....	13
8.4 Understanding the Results.....	13
8.5 Settings.....	14
9. Troubleshooting.....	15
10. Specifications.....	16

1. Disclaimer

This equipment is not a medical device and should not be used for any medical purposes or any medical conditions. This product is intended only for general wellness use.


2. Caution

- This product is NOT for people with a pacemaker.

3. Unpacking

Main Unit; Quick Guide; Charging / Data Cable;
ECG Cable; Reusable ECG Electrodes;
Carry Case; Neck Strap

4. About the Main Unit

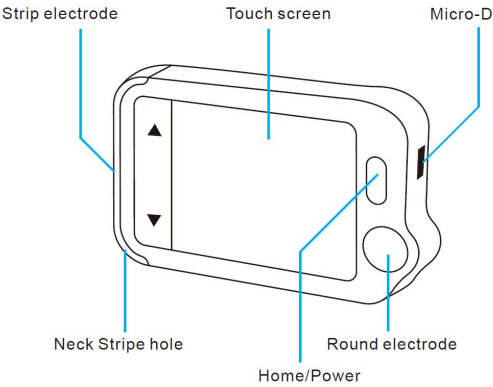5. ViHealth App

You can use **ViHealth** app to synchronize and share your data.

App name: **ViHealth**

iOS: App Store

Android: Google Play

Use your mobile phone to search for **ViHealth** in the App Store or Google Play.

Download and install the **ViHealth** mobile app.

Caution:

Uninstalling the app will delete all data!

If you need to reinstall the app, please install the new version directly, it will automatically update the old version.

6. Desktop Software

Pulsebit Browser Pro

Use the desktop software to download your Pulsebit EX data onto Windows or Mac.

Download this software from our support page:

www.getwellue.com

Menu: Support->PC Software

7. Set Up

7.1 Charging

Connect the device to a standard USB charging adapter with the provided charging cable.

7.2 Powering On/Off

Press and hold the home button for 1 second to power on, and 3 seconds to power off.

It will **turn off automatically** if it is not measuring or operating.

7.3 Initial Settings

When the device is powered on for the first time, please set the current date and time.

7.4 First Check

- 1) On the homepage, touch the **Heart Check** icon, then choose an **ECG Type Lead II (Read notes)**.

- 2) Hold the device as shown in the device display.
- 3) Remain still and relaxed while taking measurements until complete. Next, the results will be displayed.

Note:

- Moisten skin (the area that the electrode touches) *may help to get a better waveform*.
- *You can get different signal amplitudes from Lead I, Lead II, or Chest Lead. Usually, the signals from Lead II and Chest Lead are stronger*

- *For Lead I, rest your hands on a table or thigh to help you keep relax throughout the check. If your hand is trembling, try your left wrist or bicep instead of your left hand.*
- *Touching the rectangular electrode on the back of the device has little influence.*

7.5 View Results

There are 6 kinds of measurement results, Touch the " **i** " icon to view the specific definitions

Touch the "▶" icon to replay the ECG waveforms.

Touch the "C" icon to take new measurements.

You can use the **home button** at any time to return to the homepage or cancel taking measurements.

7.6 Reviewing Records

On the homepage, touch the **History** icon to review your history results.

You can choose a record to review the results and touch "▶"

icon to replay the recorded ECG waveforms.

Toch"  " icon to delete the current record.

7.7 Syncing Data to ViHealth App

Search for ViHealth in the App Store or Google Play.

Download and install the ViHealth mobile app.


- 1) On the homepage, choose **Settings->Bluetooth** to turn on the device's Bluetooth.


- 2) Switch on Bluetooth on your mobile phone.
- 3) Run ViHealth on your mobile phone. (Please **Always Allow** location and storage permissions for the Bluetooth

connection and data storage.)

- 4) Add a new device: The app will automatically search for devices and display them on the screen. Choose Pulsebit EX and wait for the data to download.

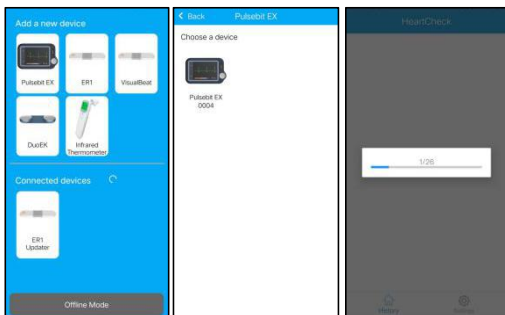- 5) Review report: Choose a single record, and you can enter a detailed report.

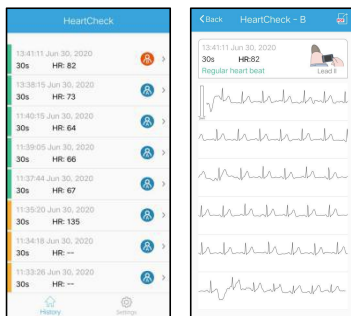
- 6) Share report: On the report details page, tap the share button in the upper right corner to share your report.

- 7) Delete a record: Choose a single record, swipe left to delete it.

Note: This operation will only delete the records in the app,

not the records on Pulsebit EX device.

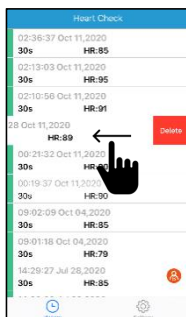7.8 Using the Software “Pulsebit Browser Pro”

Download the latest version from:

www.getwellue.com

Menu: [Support->PC Software](#)

Note:

- 1) *If you have downloaded the software before, please update it to the latest version.*
- 2) *Occasionally, firewall or anti-virus software may block “Pulsebit Browser”. Allow it.*
- 3) *Please check the help document in the software for more information.*

Download and Install *Pulsebit Browser Pro* on your computer.

To export data:

- 1) Connect the device to the PC with the **provided cable**
- 2) On the device, enter **Settings->To PC**
- 3) Run the software, click the **Download Button** to export the data on the device

With the software, you can view and print reports, which can also be exported as PDF files.

8. Exploring More Features


8.1 Choosing Different ECG Types

Before taking measurements, you can choose different ECG types (Leads). Any results will be marked with the corresponding type.


8.2 Using the ECG cable for Heart Checks

- 1) Touch <HeartCheck> to start taking measurements; **touch the icon at the bottom to switch to *External* mode.**

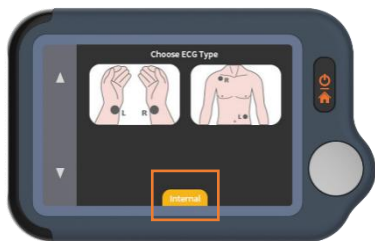- 2) Insert the ECG cable into the Micro D port, adhere the reusable electrodes onto your body as shown below.

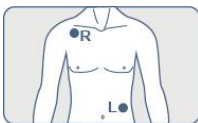★Tips:

- Use Lead II if signal is too low in Lead I
- Moisten skin (area that touch electrodes) if too dry
- Rest hands on table or thigh and keep relax
- Use the cable when you have difficulty to get stable ECG without cable

8.3 About Reusable Electrodes

Please cover the white protective paperback and seal the electrode pads in the plastic bag after each use.

We recommend replacing the electrodes after 50 uses.

Electrode performance will decrease after excessive use.

You can purchase more electrodes on Amazon or our website getwellue.com

8.4 Understanding the Results

Heart Rate	The number of times your heart beats per minute calculated from the ECG. The reference range for adults is 50-100 beats per minute in a relaxed state and higher when stimulated. Control your maximum heart rate during fitness or exercise to avoid overstressing the heart.
ECG	An electrocardiogram (ECG) is a test that records the heart’s electrical activity. Regular heart beat - the recorded ECG shows a

	usual and regular pattern. Irregular heart beat - the recorded ECG shows an erratic pattern. It may be caused by oversteering the heart or other conditions.
Poor signal	Low signal amplitude or noise. Changing to Mode II may help.
Unable to analyze	ECG recording is incomplete.

8.5 Settings

User Mode: Settings->User Mode->touch to choose

One-User/Two-User modes

When you switch to Dual User mode, the homepage will show User A and User B icons. Choose a user to start taking measurements, and the results will be marked A or B.

In single user mode, the results are not marked.

ECG Length: Settings->ECG Length, for 30s/60s/5 minutes

Note: only measurements through cables are supported in 5-minute mode

Brightness: Settings-> Touch the brightness meter to adjust the brightness

Voice: Settings->Touch the Voice switch to turn the heartbeat voice on and off.

Other settings: Date& Time, Erase Data, Factory Reset, About.

9. Troubleshooting

Problem	Possible Cause	Solution
The device does not turn on.	1. The battery may be low. 2. The device might be damaged	1. Charge the battery and try again. 2. Please contact your local distributor.
Low battery indicator is blinking	The battery is low.	Charge the battery and try again.
The ECG waveform amplitude is small	The Lead you have chosen is not suitable.	Swap with another Lead and try again.
ECG waveform drifts or disappears and “Check	1. The pressure exerted on the electrode is not stable or too great.	1. Hold the device so that it is stable. 2. Try to keep

Contact" displays.	2. Hand or body may be moving.	perfectly still and retest.
--------------------	--------------------------------	-----------------------------

10. Specifications

Size	88×56×13 mm
Weight	64 g (main unit)
Display	2.4" touch screen, color, backlight
Battery type	Rechargeable lithium-polymer battery
Battery run time	>500 checks
Charge time	2 -3 hours
Heart rate range	30 to 250 bpm
Storage	30s x 100 pieces
User	One or Two
ECG length	30s, 60s, 5 minutes

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Shenzhen Viatom Technology Co., Ltd.